

Iodine-123 Radiochemical Sodium Iodide Solution

Sodium Iodide (I-123)

Half-Life: 13.2 hours

IODINE-123 RADIOCHEMICAL SODIUM IODIDE SOLUTION

Iodine-123 is an ideal research tool for Neurology, Oncology and Cardiology.

Researchers prefer using iodine-123 because of the short half-life, clear SPECT images, and because iodine-123 is easy to handle, shield and dispose of.

Nordion has an unparalleled base of technical expertise and production capabilities. State-of-the-art facilities provide a reliable supply of isotopes to customers around the world.


PROCESS DATA

Nuclear Reaction

Primary



Secondary

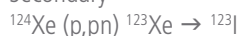

Assay

High resolution gamma spectrometry at 159 keV

PRODUCT SPECIFICATION

Chemical Form

Sodium iodide in 0.1N NaOH

Specific Activity

No carrier added

Activity Concentration*

At customer request

Radiochemical Purity

≥ 95% iodide

Radiopurity*

≥ 99.8% ^{123}I

≤ 0.2% total other nuclides

*Activity and radiopurity specifications are at customer's calibration time.

Please note that the radiochemicals listed are not tested for sterility or pyrogenicity. Verification of their suitability for use in humans is the sole responsibility of the purchaser.

PACKAGING

- Kontes vial – crimp cap glass V-vial (1 mL, 3 mL, 5 mL)
- Hypovial – crimp cap glass vial with septum (6 mL, 10 mL, 20 mL, 30 mL) (Other vial options are available).

AVAILABILITY

Dispensed Monday, Tuesday, Wednesday, Thursday

PRODUCT ORIGIN

Vancouver, Canada

WORLD CLASS DELIVERY SERVICE

Nordion has a world class distribution system that includes a unique, software-based order management system. Carefully managed transportation logistics ensure products reach their destination on time.



Iodine-123 Decay Table Half-Life: 13.2 hours

HOURS	0	0.25	0.5	0.75	1	1.25	1.5	1.75	2	2.25	2.5	2.75
0	1.0000	0.9870	0.9741	0.9614	0.9488	0.9365	0.9243	0.9122	0.9003	0.8886	0.8770	0.8655
3	0.8542	0.8431	0.8321	0.8213	0.8105	0.8000	0.7895	0.7792	0.7691	0.7591	0.7492	0.7394
6	0.7297	0.7202	0.7108	0.7016	0.6924	0.6834	0.6745	0.6657	0.6570	0.6484	0.6400	0.6316
9	0.6234	0.6152	0.6072	0.5993	0.5915	0.5838	0.5762	0.5686	0.5612	0.5539	0.5467	0.5396
12	0.5325	0.5256	0.5187	0.5120	0.5053	0.4987	0.4922	0.4858	0.4794	0.4732	0.4670	0.4609
15	0.4549	0.4490	0.4431	0.4373	0.4316	0.4260	0.4204	0.4150	0.4096	0.4042	0.3989	0.3937
18	0.3886	0.3835	0.3785	0.3736	0.3687	0.3639	0.3592	0.3545	0.3499	0.3453	0.3408	0.3363
21	0.3320	0.3276	0.3234	0.3191	0.3150	0.3109	0.3068	0.3028	0.2989	0.2950	0.2911	0.2873
24	0.2836	0.2799	0.2762	0.2726	0.2691	0.2656	0.2621	0.2587	0.2553	0.2520	0.2487	0.2454
27	0.2422	0.2391	0.2360	0.2329	0.2299	0.2269	0.2239	0.2210	0.2181	0.2152	0.2124	0.2097
30	0.2069	0.2042	0.2016	0.1989	0.1964	0.1938	0.1913	0.1888	0.1863	0.1839	0.1815	0.1791
33	0.1768	0.1745	0.1722	0.1699	0.1677	0.1655	0.1634	0.1613	0.1592	0.1571	0.1550	0.1530
36	0.1510	0.1490	0.1471	0.1452	0.1433	0.1414	0.1396	0.1378	0.1360	0.1342	0.1324	0.1307
39	0.1290	0.1273	0.1257	0.1240	0.1224	0.1208	0.1192	0.1177	0.1161	0.1146	0.1131	0.1117
42	0.1102	0.1088	0.1073	0.1059	0.1046	0.1032	0.1019	0.1005	0.0992	0.0979	0.0966	0.0954
45	0.0941	0.0929	0.0917	0.0905	0.0893	0.0882	0.0870	0.0859	0.0848	0.0836	0.0826	0.0815
48	0.0804	0.0794	0.0783	0.0773	0.0763	0.0753	0.0743	0.0734	0.0724	0.0715	0.0705	0.0696
51	0.0687	0.0678	0.0669	0.0660	0.0652	0.0643	0.0635	0.0627	0.0618	0.0610	0.0602	0.0595
54	0.0587	0.0579	0.0572	0.0564	0.0557	0.0550	0.0542	0.0535	0.0528	0.0521	0.0515	0.0508
57	0.0501	0.0495	0.0488	0.0482	0.0479	0.0469	0.0463	0.0457	0.0451	0.0445	0.0440	0.0434

The half-life values stated are as of time of publication. Nordion values may vary.

CORPORATE HEADQUARTERS

447 March Road
Ottawa, ON, Canada K2K 1X8
Tel: +1 613 592 2790
Fax: +1 613 592 6937

REGIONAL OFFICE

4004 Wesbrook Mall
Vancouver, BC, Canada V6T 2A3
Tel: +1 604 228 8952
Fax: +1 604 228 5990

ASIA PACIFIC SALES OFFICES

Hong Kong

901 Matheson Centre
3 Matheson Street
Causeway Bay, Hong Kong
Tel: +852 2827 8666
Fax: +852 2827 8302

Japan

Room 905, Tokyo Royal Plaza
1-18-11, Uchi-kanda, Chiyoda-ku
Tokyo 101-0047, Japan
Tel: +81 3 5283 6872
Fax: +81 3 5283 6873

EUROPEAN OFFICE

Zoning Industriel
Avenue de l'Espérance
B-6220 Fleurus, Belgium
Tel: +32 71 82 35 86
Fax: +32 71 82 36 66

www.nordion.com