



## Instruction Manual

# Zeocin™

Catalog nos. R250-01, R250-05

**Version K**  
October 19, 2010  
25-0078



# Table of Contents

|  |           |
|--|-----------|
| Table of Contents .....                    | iii       |
| Important Information .....                | v         |
| <b>Methods.....</b>                        | <b>1</b>  |
| Overview .....                             | 1         |
| Zeocin™ Selection in <i>E. coli</i> .....  | 4         |
| Zeocin™ Selection in Yeast.....            | 6         |
| Zeocin™ Selection in Mammalian Cells.....  | 8         |
| <b>Appendix .....</b>                      | <b>12</b> |
| HEK 293 Cells Under Zeocin™ Selection..... | 12        |
| COS Cells Under Zeocin™ Selection .....    | 13        |
| Technical Service.....                     | 14        |
| References.....                            | 16        |



# Important Information

---

## Shipping and Storage

Zeocin™ is shipped on blue ice, and is supplied as a 100 mg/ml solution in deionized, autoclaved water. **Store at -20°C.**

| Catalog No. | Amount | How Supplied |
|-------------|--------|--------------|
| R250-01     | 1 g    | 8 x 1.25 ml  |
| R250-05     | 5 g    | 50 ml        |

---

## Product Qualification

Zeocin™ is lot qualified by demonstrating that LB media containing 25 µg/ml Zeocin™ prevents growth of the *E. coli* strain, TOP10.

---



# Methods

## Overview

---


### Introduction

Zeocin™ is a member of the bleomycin / phleomycin family of antibiotics isolated from *Streptomyces*. It shows strong toxicity against bacteria, fungi (including yeast), plants and mammalian cell lines (Calmels *et al.*, 1991; Drocourt *et al.*, 1990; Gatignol *et al.*, 1987; Mulsant *et al.*, 1988; Perez *et al.*, 1989). Since Zeocin™ is active in both bacteria and mammalian cell lines, vectors can be designed that carry only one drug resistance marker for selection.

---

### Description

Zeocin™ is a formulation of phleomycin D1, a basic, water-soluble, copper-chelated glycopeptide isolated from *Streptomyces verticillus*. The presence of copper gives the solution its blue color. This copper-chelated form is inactive. When the antibiotic enters the cell, the copper cation is reduced from  $\text{Cu}^{2+}$  to  $\text{Cu}^{1+}$  and removed by sulfhydryl compounds in the cell. Upon removal of the copper, Zeocin™ is activated and will bind DNA and cleave it, causing cell death. The structure of Zeocin™ is shown below (Berdy, 1980).

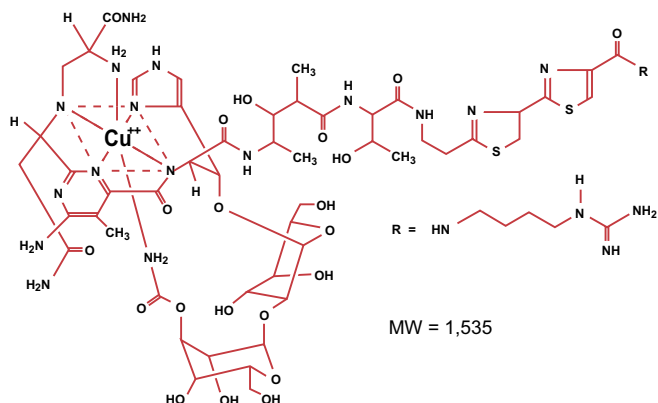

*continued on next page*

# Overview, continued

---

## Resistance to Zeocin™

A Zeocin™ resistance protein has been isolated and characterized (Calmels *et al.*, 1991; Drocourt *et al.*, 1990). This 13,665 Da protein, the product of the *Sh ble* gene (*Streptoalloteichus hindustanus* bleomycin gene), binds stoichiometrically to Zeocin™ and inhibits its DNA strand cleavage activity. Expression of this protein in eukaryotic and prokaryotic hosts confers resistance to Zeocin™.

---

## Handling Zeocin™

- High ionic strength and acidity or basicity inhibit the activity of Zeocin™. Therefore, we recommend that you reduce the salt in bacterial medium and adjust the pH to 7.5 to keep the drug active (see **Low Salt LB Medium**, page 4).
  - Store Zeocin™ at -20°C and thaw on ice before use.
  - Zeocin™ is light sensitive. Store the drug and plates or medium containing the drug in the dark.
  - Wear gloves, a laboratory coat, and safety glasses when handling Zeocin™ containing solutions.
  - Do not ingest or inhale solutions containing the drug.
  - Be sure to bandage any cuts on your fingers to avoid exposure to the drug.
- 

*continued on next page*



# Overview, continued

---

## Concentrations of Zeocin™ to Use for Selection

Zeocin™ and the *Sh ble* gene can be used for selection in mammalian cells (Drocourt *et al.*, 1990; Mulsant *et al.*, 1988), plants (Perez *et al.*, 1989), yeast (Calmels *et al.*, 1991; Gatignol *et al.*, 1987), and prokaryotes (Drocourt *et al.*, 1990). Suggested concentrations of Zeocin™ to use for selection in mammalian tissue culture cells, yeast, and *E. coli* are listed below.

| Organism        | Zeocin™ Concentration and Selective Medium |
|-----------------|--|
| <i>E. coli</i>  | 25-50 µg/ml in Low Salt LB medium*         |
| Yeast           | 50-300 µg/ml in YPD or minimal medium      |
| Mammalian cells | 50-1000 µg/ml (varies with cell line)      |

\*For efficient selection, the concentration of NaCl should not exceed 5 g/liter.

---

# Zeocin™ Selection in *E. coli*

---

## Introduction

Use 25-50 µg/ml of Zeocin™ for selection in *E. coli*. High salt and extremes in pH will inhibit the activity of Zeocin™ (see recommendations below).

---

## *E. coli* Host

Any *E. coli* strain that contains the complete Tn5 transposable element (*i.e.* DH5αF'IQ, SURE, SURE2) encodes the *ble* (bleomycin) resistance gene. These strains will confer resistance to Zeocin™. For the most efficient selection, use an *E. coli* strain that does not contain the Tn5 gene (*i.e.* TOP10, DH5, DH10, etc.).

---

## Ionic Strength and pH

Extremes in pH and high ionic strength will inhibit the activity of Zeocin™. To optimize selection in *E. coli*, the salt concentration must be < 110 mM and the pH must be 7.5. A recipe for Low Salt LB is provided below to optimize selection in *E. coli*.

---

## Low Salt LB Medium

10 g Tryptone  
5 g NaCl  
5 g Yeast Extract

1. Combine the dry reagents above and add deionized, distilled water to 950 ml. Adjust pH to 7.5 with 1 N NaOH and bring the volume up to 1 liter. For plates, add 15 g/L agar before autoclaving.
2. Autoclave on liquid cycle at 15 psi and 121°C for 20 minutes.
3. Thaw Zeocin™ on ice **and vortex before removing an aliquot.**
4. Allow the medium to cool to at least 55°C before adding the Zeocin™ to 25 µg/ml final concentration.

Store plates and unused medium at +4°C in the dark. Plates and medium containing Zeocin™ are stable for 1-2 weeks.

---

*continued on next page*

# Zeocin™ Selection in *E. coli*, continued

---



imMedia™ is available from Invitrogen for fast and easy microwaveable preparation of Low Salt LB medium or agar containing Zeocin™. See below for ordering information. For more information, see our Web site ([www.invitrogen.com](http://www.invitrogen.com)) or call Technical Service (see page 14).

| Medium              | Quantity    | Catalog no. |
|---------------------|-------------|-------------|
| imMedia™ Zeo Liquid | 20 pouches† | Q620-20     |
| imMedia™ Zeo Agar   | 20 pouches† | Q621-20     |

†Each pouch provides sufficient reagents to prepare 200 ml of liquid medium or 8-10 standard size agar plates.

---

# Zeocin™ Selection in Yeast

---

## Introduction

We have successfully transformed plasmids conferring Zeocin™ resistance into *Saccharomyces cerevisiae* and *Pichia pastoris*. The concentration of Zeocin™ required to select resistant transformants may range from 50 to 300 µg/ml, depending on the strain, pH, and ionic strength. Guidelines are provided below to assist you with selecting Zeocin™-resistant transformants.

---



## Important

We do not recommend spheroplasting for transformation of yeast with plasmids containing the Zeocin™ resistance marker. Spheroplasting involves removal of the cell wall to allow DNA to enter the cell. Cells must first regenerate the cell wall before they are able to express the Zeocin™ resistance gene. Plating spheroplasts directly onto selective medium containing Zeocin™ will result in complete cell death.

---

## Transformation Method

We recommend electroporation, lithium cation protocols, or our EasyComp™ Kits for transformation of yeast with vectors that encode resistance to Zeocin™. Electroporation yields 10<sup>3</sup> to 10<sup>4</sup> transformants per µg of linearized DNA and does not destroy the cell wall of yeast. If you do not have access to an electroporation device, use chemical methods or one of the EasyComp™ Kits listed below.

| Kit  | Reactions              | Catalog no. |
|--|------------------------|-------------|
| <i>S. c.</i> EasyComp™ Transformation Kit (for <i>Saccharomyces cerevisiae</i> ) | 6 x 20 transformations | K5050-01    |
| <i>Pichia</i> EasyComp™ Transformation Kit (for <i>Pichia pastoris</i> )         | 6 x 20 transformations | K1730-01    |

*continued on next page*

# Zeocin™ Selection in Yeast, continued

---

## Ionic Strength and pH

Since yeast vary in their susceptibility to Zeocin™, we recommend that you perform a kill curve to determine the lowest concentration of Zeocin™ needed to kill the untransformed host strain. In addition, the pH of the selection medium may affect the concentration of Zeocin™ needed to select resistant transformants. We recommend that you test media adjusted to different pH values (6.5 to 8) for the one that allows you to use the lowest possible concentration of Zeocin™.

---

## Selection in Yeast

For successful selection of Zeocin™-resistant transformants we recommend the following:

- After transformation (either by electroporation or chemical transformation), allow the cells to recover for one hour in YPD medium.
  - For electroporated cells, plate your transformants on YPD containing 1 M sorbitol. Sorbitol allows better recovery of the cells after electroporation.
  - For chemically transformed cells, plate cells on YPD or minimal plates.
  - Plate several different volumes (*i.e.* 10, 25, 50, 100, and 200 µl) of the transformation reaction. Plating at low cell densities favors efficient Zeocin™ selection.
-

# Zeocin™ Selection in Mammalian Cells

---

## Introduction

Mammalian cells exhibit a wide range of susceptibility to Zeocin™. Concentrations of Zeocin™ used to select stable cell lines may range from 50 to 1000 µg/ml, with the average being around 250 to 400 µg/ml. Factors that affect selection include ionic strength, cell line, cell density, and growth rate. Review the guidelines below to ensure successful selection of stable cell lines.

---



## Important

Zeocin™'s mechanism of killing is quite different from neomycin and hygromycin. **Cells do not round up and detach from the plate.** Sensitive cells may exhibit the following morphological changes upon exposure to Zeocin™:

- Vast increase in size (similar to the effects of cytomegalovirus infecting permissive cells)
- Abnormal cell shape, including the appearance of long appendages
- Presence of large empty vesicles in the cytoplasm (breakdown of the endoplasmic reticulum and Golgi apparatus or scaffolding proteins)
- Breakdown of plasma and nuclear membrane (appearance of many holes in these membranes)

Eventually, these cells will completely break down and only cellular debris will remain.

Zeocin™-resistant cells should continue to divide at regular intervals to form distinct colonies. There should not be any distinct morphological changes in Zeocin™-resistant cells when compared to cells not under selection with Zeocin™.

---

## Examples

To see photographs of HEK 293 and COS1 cells undergoing selection in the presence of Zeocin™, refer to the **Appendix**, pages 12 and 13.

---

*continued on next page*

# Zeocin™ Selection in Mammalian Cells, continued

---

## Ionic Strength and pH

For selection in mammalian cells, physiological ionic strength and pH are much more important for cell growth, so more Zeocin™ may be needed for selection relative to yeast or bacteria.

---

## Selection in Mammalian Cell Lines

To generate a stable cell line that expresses your protein from an expression construct, you need to determine the minimum concentration required to kill your untransfected host cell line (see **Determination of Zeocin™ Sensitivity**, below). In general, it takes 2-6 weeks to generate foci with Zeocin™, depending on the cell line. Because individual cells can express protein at varying levels, it is important to isolate several foci to expand into stable cell lines.

---

## Determining Zeocin™ Sensitivity

Determine the minimal concentration of Zeocin™ required to kill the untransfected parental cell line using the protocol below.

1. Plate or split a confluent plate so the cells will be approximately 25% confluent. Prepare a set of 8 plates. Grow cells for 24 hours.
2. Remove medium and then add medium with varying concentrations of Zeocin™ (0, 50, 100, 200, 400, 600, 800, and 1000 µg/ml) to each plate.
3. Replenish the selective medium every 3-4 days and observe the percentage of surviving cells over time. Select the concentration that kills the majority of the cells in the desired number of days (within 1-2 weeks).

If you have trouble distinguishing viable cells by observation, we recommend counting the number of viable cells by trypan blue exclusion to determine the appropriate concentration of Zeocin™ required to prevent growth.

---

*continued on next page*

# Zeocin™ Selection in Mammalian Cells, continued

---

## Selection Tip

Some cells may be more resistant to Zeocin™ than others. If cells are rapidly dividing, Zeocin™ may not be effective at low concentrations. We suggest trying the following protocol to overcome this resistance:

1. Split cells into medium containing Zeocin™.
2. Incubate cells at 37°C for 2-3 hours until the cells have attached to the culture dish.
3. Remove the plates from the incubator and place the cells at +4°C for 2 hours. Be sure to buffer the medium with HEPES.
4. Return the cells to 37°C.

Incubating the cells at +4°C will stop the cell division process for a short time, allow Zeocin™ to act, and result in cell death.

---

*continued on next page*



# Zeocin™ Selection in Mammalian Cells, continued

---

## Selecting Stable Integrants

Once you have determined the appropriate Zeocin™ concentration to use for selection, you can generate a stable cell line with your construct.

1. Transfect your cell line and plate onto 100 mm culture plates. Include a sample of untransfected cells as a negative control.
2. After transfection, wash the cells once with 1X PBS and add fresh medium to the cells.
3. Forty-eight to 72 hours after transfection, split the cells using various dilutions into fresh medium containing Zeocin™ at the pre-determined concentration required for your cell line. By using different dilutions, you will have a better chance at identifying and selecting foci.

**Note:** If your cells are more resistant to Zeocin™, you may want to use the selection tip described on the previous page. Simply split cells into medium containing Zeocin™, incubate at 37°C for 2-3 hours to let cells attach, then place the cells at +4°C for 2 hours. Remember to buffer the medium with HEPES.

4. Feed the cells with selective medium every 3-4 days until cell foci are identified.
  5. Pick and transfer colonies to either 96- or 48-well plates. Grow cells to near confluence before expanding to larger wells or plates.
- 

## Maintaining Stable Cell Lines

To maintain stable cell lines, you may:

- Maintain the cells in the same concentration of Zeocin™ you used for selection
  - Reduce the concentration of Zeocin™ by half
  - Reduce the concentration of Zeocin™ to the concentration that just prevents growth of sensitive cells but does not kill them (refer to your kill curve experiment)
-

# Appendix

## HEK 293 Cells Under Zeocin™ Selection

---

### Introduction

The photographs below show HEK 293 cells (Graham *et al.*, 1977) undergoing Zeocin™ selection. Cells were cultured in DMEM containing 10% FBS, 1 mM L-glutamine, and 400 µg/ml Penicillin-Streptomycin in the absence or presence of 400 µg/ml Zeocin™.

**Panel A:** 293 cells not exposed to Zeocin™

**Panels B and C:** 293 cells after 3 days in selective medium




**A. Unselected cells**



**B. Zeocin™ -sensitive cells**

Long appendages may appear to grow out from the cell as the plasma membrane breaks down (see filled arrow in this panel and below).

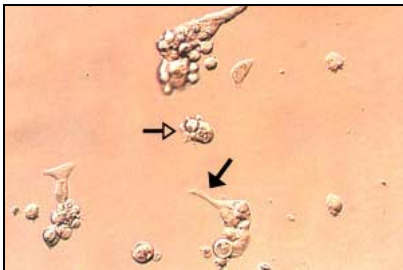

**C. Zeocin™ -sensitive cells**

Cells will begin to disintegrate and cell particles may be observed in the medium (see open arrow in this panel).

---

# COS Cells Under Zeocin™ Selection

---

## Introduction

The photographs below show COS1 cells undergoing Zeocin™ selection. Cells were cultured in DMEM containing 10% FBS, 1 mM L-glutamine, and 400 µg/ml Penicillin-Streptomycin in the absence or presence of 400 µg/ml Zeocin™.

**Panel A:** COS1 cells not exposed to Zeocin™

**Panels B and C:** COS1 cells after 3 days in selective medium



**Unselected cells**



**Zeocin™-sensitive cells**

Cells will begin to disintegrate and cell particles may be observed in the medium (see filled arrow in this panel).



**Zeocin™-sensitive cells**

Long appendages may appear to grow out from the cell as the plasma membrane breaks down (see filled arrow in this panel).

---

# Technical Service

---

## World Wide Web



Visit the Invitrogen Web Resource using your World Wide Web browser. At the site, you can:

- Get the scoop on our hot new products and special product offers
- View and download vector maps and sequences
- Download manuals in Adobe® Acrobat® (PDF) format
- Explore our catalog with full color graphics
- Obtain citations for Invitrogen products
- Request catalog and product literature

Once connected to the Internet, launch your web browser (Internet Explorer 5.0 or newer or Netscape 4.0 or newer), then enter the following location (or URL):

**<http://www.invitrogen.com>**

...and the program will connect directly. Click on underlined text or outlined graphics to explore. Don't forget to put a bookmark at our site for easy reference!

---

## Contact Us

For more information or technical assistance, call, write, fax, or email. Additional international offices are listed on our web page ([www.invitrogen.com](http://www.invitrogen.com)).

### Corporate Headquarters:

Invitrogen Corporation  
1600 Faraday Avenue  
Carlsbad, CA 92008 USA  
Tel: 1 760 603 7200  
Tel (Toll Free): 1 800 955 6288  
Fax: 1 760 602 6500  
E-mail:  
[tech\\_service@invitrogen.com](mailto:tech_service@invitrogen.com)

### European Headquarters:

Invitrogen Ltd  
Inchinnan Business Park  
3 Fountain Drive  
Paisley PA4 9RF, UK  
Tel: +44 (0) 141 814 6100  
Tech Fax: +44 (0) 141 814 6117  
E-mail: [eurotech@invitrogen.com](mailto:eurotech@invitrogen.com)

---

*continued on next page*

# Technical Service, continued

---

## MSDS Requests

To request an MSDS, visit our Web site ([www.invitrogen.com](http://www.invitrogen.com)).

1. On the home page, go to the left-hand column under 'Technical Resources' and select 'MSDS Requests'.
  2. Follow instructions on the page and fill out all the required fields.
  3. To request additional MSDSs, click the 'Add Another' button.
  4. All requests will be faxed unless another method is selected.
  5. When you are finished entering information, click the 'Submit' button. Your MSDS will be sent within 24 hours.
- 

## Limited Warranty

Invitrogen is committed to providing our customers with high-quality goods and services. Our goal is to ensure that every customer is 100% satisfied with our products and our service. If you should have any questions or concerns about an Invitrogen product or service, please contact our Technical Service Representatives.

Invitrogen warrants that all of its products will perform according to the specifications stated on the certificate of analysis. The company will replace, free of charge, any product that does not meet those specifications. This warranty limits Invitrogen Corporation's liability only to the cost of the product. No warranty is granted for products beyond their listed expiration date. No warranty is applicable unless all product components are stored in accordance with instructions. Invitrogen reserves the right to select the method(s) used to analyze a product unless Invitrogen agrees to a specified method in writing prior to acceptance of the order.

Invitrogen makes every effort to ensure the accuracy of its publications, but realizes that the occasional typographical or other error is inevitable. Therefore Invitrogen makes no warranty of any kind regarding the contents of any publications or documentation. If you discover an error in any of our publications, please report it to our Technical Service Representatives.

**Invitrogen assumes no responsibility or liability for any special, incidental, indirect or consequential loss or damage whatsoever. The above limited warranty is sole and exclusive. No other warranty is made, whether expressed or implied, including any warranty of merchantability or fitness for a particular purpose.**

---

# References

---

- Berdy, J. (1980) Bleomycin-Type Antibiotics. In *Amino Acid and Peptide Antibiotics*, J. Berdy, ed. (Boca Raton, FL: CRC Press), pp. 459-497.
- Calmels, T., Parriche, M., Burand, H., and Tiraby, G. (1991). High Efficiency Transformation of *Tolypocladium geodes* Conidiospores to Phleomycin Resistance. *Curr. Genet.* 20, 309-314.
- Drocourt, D., Calmels, T. P. G., Reynes, J. P., Baron, M., and Tiraby, G. (1990). Cassettes of the *Streptoalloteichus hindustanus ble* Gene for Transformation of Lower and Higher Eukaryotes to Phleomycin Resistance. *Nucleic Acids Res.* 18, 4009.
- Gatignol, A., Baron, M., and Tiraby, G. (1987). Phleomycin Resistance Encoded by the *ble* Gene from Transposon Tn5 as a Dominant Selectable Marker in *Saccharomyces cerevisiae*. *Molecular and General Genetics* 207, 342-348.
- Graham, F. L., Smiley, J., Russell, W. C., and Nairn, R. (1977). Characteristics of a Human Cell Line Transformed by DNA from Human Adenovirus Type 5. *J. Gen. Virol.* 36, 59-74.
- Mulsant, P., Tiraby, G., Kallerhoff, J., and Perret, J. (1988). Phleomycin Resistance as a Dominant Selectable Marker in CHO Cells. *Somat. Cell Mol. Genet.* 14, 243-252.
- Perez, P., Tiraby, G., Kallerhoff, J., and Perret, J. (1989). Phleomycin Resistance as a Dominant Selectable Marker for Plant Cell Transformation. *Plant Mol. Biol.* 13, 365-373.

---

©1998-2002, 2010 Invitrogen Corporation. All rights reserved.





**Corporate Headquarters:**

*Invitrogen Corporation  
1600 Faraday Avenue  
Carlsbad, California 92008  
Tel: 1 760 603 7200  
Tel (Toll Free): 1 800 955 6288  
Fax: 1 760 603 7229  
Email: tech\_service@invitrogen.com*

**European Headquarters:**

*Invitrogen Ltd  
3 Fountain Drive  
Inchinnan Business Park  
Paisley PA4 9RF, UK  
Tel (Free Phone Orders): 0800 269 210  
Tel (General Enquiries): 0800 5345 5345  
Fax: +44 (0) 141 814 6287  
Email: eurotech@invitrogen.com*

**International Offices:**

*Argentina 5411 4556 0844  
Australia 1 800 331 627  
Austria 0800 20 1087  
Belgium 0800 14894  
Brazil 0800 11 0575  
Canada 800 263 6236  
China 10 6849 2578  
Denmark 80 30 17 40*

*France 0800 23 20 79  
Germany 0800 083 0902  
Hong Kong 2407 8450  
India 11 577 3282  
Italy 02 98 22 201  
Japan 03 3663 7974  
The Netherlands 0800 099 3310  
New Zealand 0800 600 200  
Norway 00800 5456 5456*

*Spain & Portugal 900 181 461  
Sweden 020 26 34 52  
Switzerland 0800 848 800  
Taiwan 2 2651 6156  
UK 0800 838 380  
For other countries see our Web site*

[www.invitrogen.com](http://www.invitrogen.com)

