

Daily QATM 3

Daily QA Accelerated

The most advanced and easy to use instrument for routine quality assurance of radiotherapy beams.

U.S. Patent No. 6,125,335



rf-Daily QA 3
(wireless)

Daily QA 3

Benefits

- > Utilizes both ion chambers and diodes for optimal results
- > Simultaneously check:
 - Output, flatness, symmetry, field size, energy
- > Deliver rotational beams if desired
- > FFF beams supported
- > No flipping or additional buildup required
- > Trend analysis database included
- > Real-time wireless version available (rf-Daily QA 3):
 - Keeps cables off treatment room floor for enhanced safety

Hardware Features


- > 13 ion chambers:
 - Output, flatness, symmetry, energy
- > 12 SunPoint[®] Diode Detectors:
 - Measure light-radiation field coincidence
- > Automatic temperature and pressure correction
- > Integrated buildup; no additional buildup required
- > Real-time measurements

Software Features


- > Simple two-step operation - 'Start' and then 'Record'
- > Real-time display of measured data
- > Green, amber or red light indicates pass, warning, or fail
- > Sophisticated analysis options with the click of a mouse
- > Use different devices for a template without creating a new baseline
- > Exports PDF reports

Optional Accessory


- > Use the IMF to position the Daily QA 3 at isocenter and measure at any gantry angle

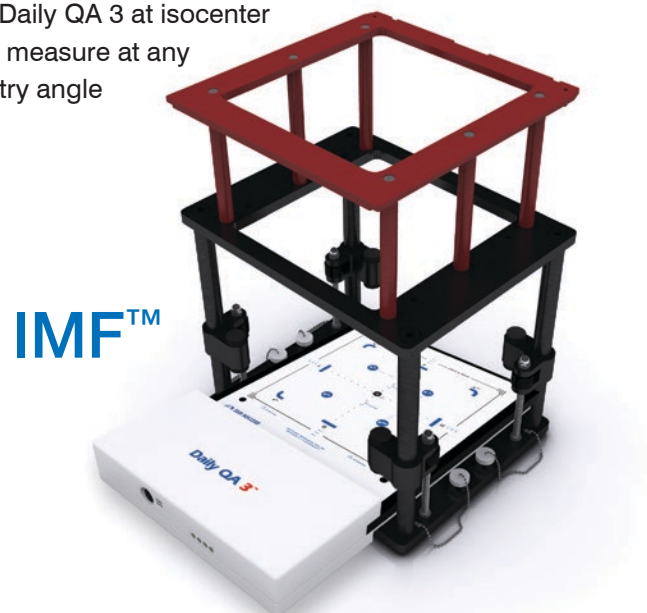

IMFTM

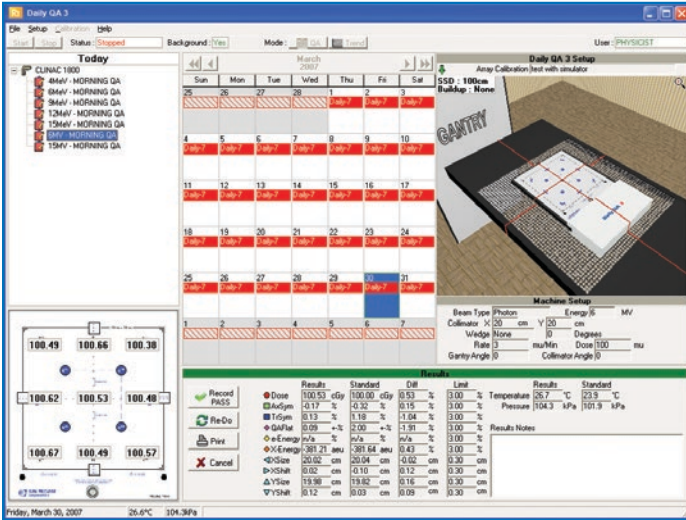


After a measurement is accepted, all collected data is saved for statistical analysis and reporting. Recording data will also initiate the next test. Users can conduct a new test, scheduled or unscheduled, at any time and data is saved in the database.

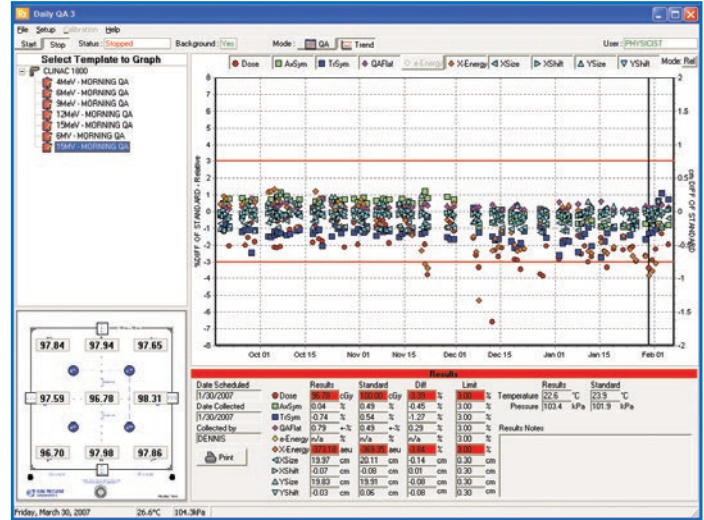

Daily QA 3 software provides a graphical presentation of data for each template. Users have the option to choose which data parameters to view and over what period to view them. A single measurement instance can be examined by clicking on a data point which corresponds to a particular day.

Features

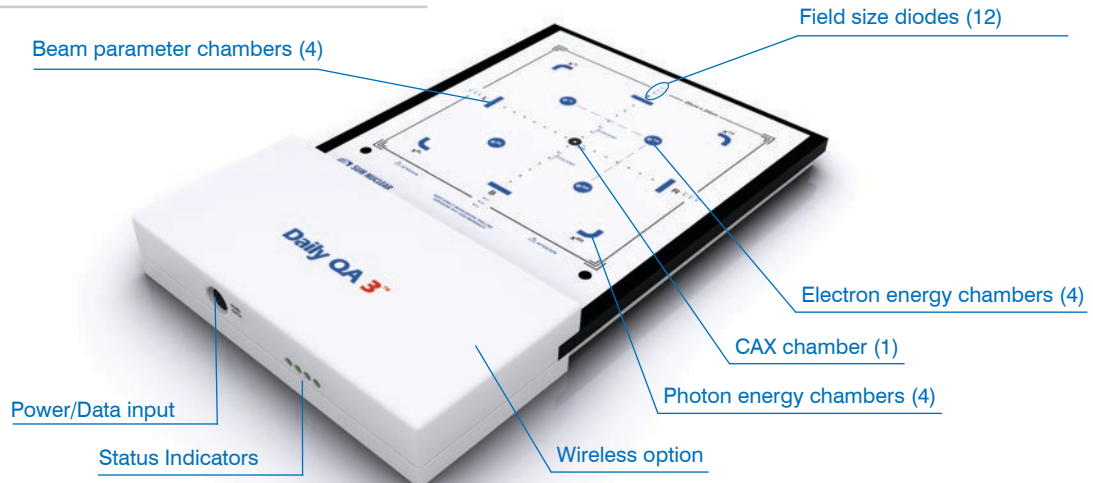

Specifications

Detector type:	SunPoint® Diode Detectors Vented Ion Chambers	Inherent backscatter (g/cm ²):	2.3
Detector quantity:	Diodes: 12 total Chambers: 13 total 1 CAX	Electron energy attenuation:	Air, Cu, Al, Fe
Detector spacing (mm):	Diodes: 5.0	Radiation measured:	• Electrons, 4MeV to 25MeV • Photons, Co-60 to 25MV
Chamber active area (cm ²):	Electron: 0.6 Photon: 0.3	rf Frequency (rf-Daily QA 3):	2.400 to 2.485GHz
Field size (cm):	20.0 x 20.0	Operating system:	Windows 2000, XP, Vista, 7 (32-64bit)
Inherent buildup (g/cm ²):	Chambers: 1.0 ± 0.1	Dimensions / Weight:	25.6 x 40.8 x 4.6cm / 5.7kg
		Number of connection cables:	Single power / data cable

