

See discussions, stats, and author profiles for this publication at: <https://www.researchgate.net/publication/346483715>

External peer review of the RTPCR test to detect SARS-CoV-2 reveals 10 major scientific flaws at the molecular and methodological level: consequences for false positive results

Preprint · November 2020

DOI: 10.5281/zenodo.4298004

CITATIONS

0

READS

66,812

22 authors, including:



Pieter Berger

W & W Research

83 PUBLICATIONS 2,010 CITATIONS

SEE PROFILE



Klaus Steger

Justus-Liebig-Universität Gießen

247 PUBLICATIONS 7,842 CITATIONS

SEE PROFILE



Stefano Scoglio

Università degli Studi di Urbino "Carlo Bo"

23 PUBLICATIONS 355 CITATIONS

SEE PROFILE



Rainer J Klement

Leopoldina Hospital, Schweinfurt

131 PUBLICATIONS 1,680 CITATIONS

SEE PROFILE

Some of the authors of this publication are also working on these related projects:



Science against panic in the COVID-19 crisis [View project](#)



THE PALEO PROJECT [View project](#)

External peer review of the RTPCR test to detect SARS-CoV-2 reveals 10 major scientific flaws at the molecular and methodological level: consequences for false positive results.

Pieter Borger ^{1*}, Rajesh K. Malhotra ², Michael Yeadon ³, Clare Craig ⁴, Kevin McKernan ⁵
Klaus Steger ⁶, Paul McSheehy ⁷, Lidiya Angelova ⁸, Fabio Franchi ⁹, Thomas Binder ¹⁰
Henrik Ullrich ¹¹, Makoto Ohashi ¹², Stefano Scoglio ¹³, Marjolein Doesburg-van Kleffens ¹⁴
Dorothea Gilbert ¹⁵, Rainer J. Klement ¹⁶, Ruth Schruefer ¹⁷, Berber W. Pieksma ¹⁸, Jan Bonte ¹⁹,
Bruno H. Dalle Carbonara ²⁰, Kevin P. Corbett ²¹, Ulrike Kämmerer ²².

* Corresponding author

ABSTRACT

In the publication entitled “Detection of 2019 novel coronavirus (2019-nCoV) by real-time RT-PCR” (Eurosurveillance 25(8) 2020) the authors present a diagnostic workflow and RT-qPCR protocol for detection and diagnostics of 2019-nCoV (now known as SARS-CoV-2), which they claim to be validated, as well as being a *robust diagnostic methodology for use in public-health laboratory settings*.

In light of all the consequences resulting from this very publication for societies worldwide, a group of independent researchers performed a point-by-point review of the aforesaid publication in which 1) all components of the presented test design were cross checked, 2) the RT-qPCR protocol-recommendations were assessed w.r.t. good laboratory practice, and 3) parameters examined against relevant scientific literature covering the field.

The published RT-qPCR protocol for detection and diagnostics of 2019-nCoV and the manuscript suffer from numerous technical and scientific errors, including insufficient primer design, a problematic and insufficient RT-qPCR protocol, and the absence of an accurate test validation. Neither the presented test nor the manuscript itself fulfils the requirements for an acceptable scientific publication. Further, serious conflicts of interest of the authors are not mentioned. Finally, the very short timescale between submission and acceptance of the publication (24 hours) signifies that a systematic peer review process was either not performed here, or of problematic poor quality.

We provide compelling evidence of several scientific inadequacies, errors and flaws. Considering the scientific and methodological blemishes presented here, we are confident that the editorial board of Eurosurveillance has no other choice but to retract the publication.

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

CONCISE REVIEW REPORT

This paper will show numerous serious flaws in the Corman-Drosten paper, the significance of which has led to worldwide misdiagnosis of infections attributed to SARS-CoV-2 and associated with the disease COVID-19. We are confronted with stringent lockdowns which have destroyed many people's lives and livelihoods, limited access to education and these imposed restrictions by governments around the world are a direct attack on people's basic rights and their personal freedoms, resulting in collateral damage for entire economies on a global scale.

There are ten fatal problems with the Corman-Drosten paper which we will outline and explain in greater detail in the following sections.

The first and major issue is that the *novel* Coronavirus SARS-CoV-2 (in the publication named 2019-nCoV and in February 2020 named SARS-CoV-2 by an international consortium of virus experts) is based on *in silico* (theoretical) sequences, supplied by a laboratory in China [1], because at the time neither control material of infectious ("live") or inactivated SARS-CoV-2 nor isolated genomic RNA of the virus was available to the authors. To date no validation has been performed by the authorship based on isolated SARS-CoV-2 viruses or full length RNA thereof.

According to Corman et al.: "*We aimed to develop and deploy robust diagnostic methodology for use in public health laboratory settings without having virus material available.*" [1]

The focus here should be placed upon the two stated aims: a) *development* and b) *deployment* of a *diagnostic test for use in public health laboratory settings*. These aims are not achievable without having any actual virus material available (e.g. for determining the infectious viral load). In any case, only a protocol with maximal accuracy can be the mandatory and primary goal in any scenario-outcome of this magnitude. Critical viral load determination is mandatory information, and it is in Christian Drosten's group responsibility to perform these experiments and provide the crucial data.

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

Nevertheless these *in silico* sequences were used to develop a RT-PCR test methodology to identify the aforesaid virus. This model was based on the assumption that the *novel* virus is very similar to SARS-CoV from 2003 (Hereafter named SARS-CoV-1) as both are beta-coronaviruses.

The PCR test was therefore designed using the genomic sequence of SARS-CoV-1 as a control material for the Sarbeco component; we know this from our personal email-communication with [2] one of the co-authors of the Corman-Drosten paper. This method to model SARS-CoV-2 was described in the Corman-Drosten paper as follows:

“the establishment and validation of a diagnostic workflow for 2019-nCoV screening and specific confirmation, designed in absence of available virus isolates or original patient specimens. Design and validation were enabled by the close genetic relatedness to the 2003 SARS-CoV, and aided by the use of synthetic nucleic acid technology.”

In short, a design relying merely on close genetic relatives does not fulfill the aim for a “robust diagnostic test” as cross reactivity and therefore false-positive results will inevitably occur.

Validation was only done in regards to *in silico* (theoretical) sequences and within the laboratory-setting, and not as required for in-vitro diagnostics with isolated genomic viral RNA. This very fact hasn't changed even after 10 months of introduction of the test into routine diagnostics.

There are numerous other severe scientific errors regarding the biomolecular design of the primers, the PCR method, as well as the molecular validation of the PCR products and methods described in the Corman-Drosten paper which are examined in detail in the following chapters. The paper itself already signifies that a large number of false positive results are generated by this test, even under controlled laboratory conditions, making it completely unsuitable as a reliable virus screening method for entire populations in an ongoing pandemic. Given the far-reaching implications, including quarantine measures, lockdowns, curfews and impacts on education etc., this paper must be immediately retracted.

DESIGN AND ERRORS in RT-PCR

The Reverse Transcription-Polymerase Chain Reaction (RT-PCR) is an important biomolecular technology to rapidly detect rare RNA fragments, which are known in advance. In the first step, RNA molecules present in the sample are reverse transcribed to yield cDNA. The cDNA is then amplified in the polymerase chain reaction using a specific primer pair and a thermostable DNA polymerase enzyme. The technology is highly sensitive and its detection limit is theoretically 1 molecule of cDNA. The specificity of the PCR is highly influenced by biomolecular design errors.

What is important when designing an RT-PCR Test and the quantitative RT-qPCR test described in the Corman-Drosten publication?

1. The primers and probes:
 - a) the concentration of primers and probes must be of optimal range (100-200 nM)
 - b) must be specific to the target-gene you want to amplify
 - c) must have an optimal percentage of GC content relative to the total nitrogenous bases (minimum 40%, maximum 60%)
 - d) for virus diagnostics at least 3 primer pairs must detect 3 viral genes (preferably as far apart as possible in the viral genome)
2. The temperature at which all reactions take place:
 - a) DNA melting temperature (>92°)
 - b) DNA amplification temperature (TaqPol specific)
 - c) T_m; the annealing temperature (the temperature at which the primers and probes reach the target binding/detachment, not to exceed 2°C per primer pair).
T_m heavily depends on GC content of the primers
3. The number of amplification cycles (less than 35; preferably 25-30 cycles); In case of virus detection, >35 cycles only detects signals which do not correlate with infectious virus as determined by isolation in cell culture [reviewed in 2]; if someone is tested by PCR as positive when a threshold of 35 cycles or higher is used (as is the case in most laboratories in Europe & the US), the probability that said person is actually infected is less than 3%, the probability that said result is a false positive is 97%

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

[reviewed in 3]

4. Molecular biological validations; amplified PCR products must be validated either by running the products in a gel with a DNA ruler, or by direct DNA sequencing
5. Positive and negative controls should be specified to confirm/refute specific virus detection
6. There should be a Standard Operational Procedure (SOP) available, which unequivocally specifies the above parameters, so that all laboratories are able to set up the exact same test conditions. To have a validated universal SOP is essential, because it enables the comparison of data within and between countries.

MINOR CONCERNS WITH THE CORMAN-DROSTEN PAPER

1. In Table 1 of the Corman-Drosten paper, different abbreviations are stated - "nM" is specified, "nm" isn't. Further in regards to correct nomenclature, nm means "nanometer" therefore nm should read nM here.
2. It is the general consensus to write genetic sequences always in the 5'-3' direction, including the reverse primers. It is highly unusual to do alignment with reverse complementary writing of the primer sequence as the authors did in figure 2 of the Corman-Drosten paper. Here, in addition, a wobble base is marked as "y" without description of the bases the Y stands for.
3. Two misleading pitfalls in the Corman-Drosten paper are that their Table 1 does not include T_m-values (annealing-temperature values), neither does it show GC-values (number of G and C in the sequences as %-value of total bases).

MAJOR CONCERNS WITH THE CORMAN-DROSTEN PAPER

A) BACKGROUND

The authors introduce the background for their scientific work as: *“The ongoing outbreak of the recently emerged novel coronavirus (2019-nCoV) poses a challenge for public health laboratories as virus isolates are unavailable while there is growing evidence that the outbreak is more widespread than initially thought, and international spread through travelers does already occur”*.

According to BBC News [4] and Google Statistics [5] there were 6 deaths world-wide on January 21st 2020 - the day when the manuscript was submitted. Why did the authors assume a challenge for public health laboratories while there was no substantial evidence at that time to indicate that the outbreak was more widespread than initially thought?

As an aim the authors declared to develop and deploy robust diagnostic methodology for use in public health laboratory settings without having virus material available. Further, they acknowledge that *“The present study demonstrates the enormous response capacity achieved through coordination of academic and public laboratories in national and European research networks.”*

B) Methods and Results

1. Primer & Probe Design

1a) Erroneous primer concentrations

Reliable and accurate PCR-test protocols are normally designed using between 100 nM and 200 nM per primer [7]. In the Corman-Drosten paper, we observe unusually high and varying primer concentrations for several primers (table 1). For the RdRp_SARSr-F and RdRp_SARSr-R primer pairs, 600 nM and 800 nM are described, respectively. Similarly, for the N_Sarbeco_F and N_Sarbeco_R primer set, they advise 600 nM and 800 nM, respectively [1]. It should be clear that these concentrations are far too high to be optimal for specific amplifications of target genes. *There exists no specified reason to use these extremely high*

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

concentrations of primers in this protocol. Rather, these concentrations lead to increased unspecific binding and PCR product amplification.

Table 1: Primers and probes (adapted from Corman-Drosten paper; erroneous primer concentrations are highlighted)

Assay/use	Oligonucleotide	Sequence ^a	Concentration ^b
RdRP gene	RdRp_SARsR-F	GTGARATGGTCATGTGGCGG	Use 600 nM per reaction
	RdRp_SARsR-P2	FAM-CAGGTGGAACCTCATCAGGAGATGC-BBQ	Specific for 2019-nCoV, will not detect SARS-CoV. Use 100 nM per reaction and mix with P1
	RdRp_SARsR-P1	FAM-CCAGGTGGWACRTCATCMGGTGATGC-BBQ	Pan Sarbeco-Probe will detect 2019-nCoV, SARS-CoV and bat-SARS-related CoVs. Use 100 nM per reaction and mix with P2
	RdRp_SARsR-R	CARATGTAAASACACTATTAGCATA	Use 800 nM per reaction
E gene	E_Sarbeco_F	ACAGGTACGTTAATAGTTAATAGCGT	Use 400 nM per reaction
	E_Sarbeco_P1	FAM-ACACTAGCCATCCTACTGCGCTTCG-BBQ	Use 200 nm per reaction
	E_Sarbeco_R	ATATTGCAGCAGTACGCACACA	Use 400 nm per reaction
N gene	N_Sarbeco_F	CACATTGGCACCCGAATC	Use 600 nm per reaction
	N_Sarbeco_P	FAM-ACTTCCTCAAGGAACAACATTGCCA-BBQ	Use 200 nm per reaction
	N_Sarbeco_R	GAGGAACGAGAAGAGGCTTG	Use 800 nm per reaction

^a W is A/T; R is G/A; M is A/C; S is G/C. FAM: 6-carboxyfluorescein; BBQ: blackberry quencher.
^b Optimised concentrations are given in nanomol per litre (nM) based on the final reaction mix, e.g. 1.5 µL of a 10 µM primer stock solution per 25 µL total reaction volume yields a final concentration of 600 nM as indicated in the table.

1b) Unspecified (“Wobbly”) primer and probe sequences

To obtain reproducible and comparable results, it is essential to distinctively define the primer pairs. In the Corman-Drosten paper we observed six unspecified positions, indicated by the letters R, W, M and S (Table 2). The letter W means that at this position there can be either an A or a T; R signifies there can be either a G or an A; M indicates that the position may either be an A or a C; the letter S indicates there can be either a G or a C on this position.

This high number of variants not only is unusual, but it also is highly confusing for laboratories. These six unspecified positions could easily result in the design of several different alternative primer sequences which do not relate to SARS-CoV-2 (2 distinct RdRp_SARsR_F primers + 8 distinct RdRp_SARS_P1 probes + 4 distinct RdRp_SARsR_R). The design variations will inevitably lead to results that are not even SARS-CoV-2 related. Therefore, the confusing unspecific description in the Corman-Drosten paper is not suitable as a Standard Operational Protocol. These unspecified positions should have been designed unequivocally.

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

These wobbly sequences have already created a source of concern in the field and resulted in a Letter to the Editor authored by Pillonel *et al.* [8] regarding blatant errors in the described sequences. These errors are self-evident in the Corman *et al.* supplement as well.

Table 2: Primers and probes (adapted from Corman-Drosten paper; unspecified (“Wobbly”) nucleotides in the primers are highlighted)

Assay/use	Oligonucleotide	Sequence ^a	Concentration ^b
RdRP gene	RdRp_SARSr-F	GTGARATGGTCATGTGTGGCGG	Use 600 nM per reaction
	RdRp_SARSr-P2	FAM-CAGGTGGAACCTCATCAGGAGATGC-BBQ	Specific for 2019-nCoV, will not detect SARS-CoV. Use 100 nM per reaction and mix with P1
	RdRp_SARSr-P1	FAM-CCAGGTGGWACRTCATCMGGTGATGC-BBQ	Pan Sarbeco-Probe will detect 2019-nCoV, SARS-CoV and bat-SARS-related CoVs. Use 100 nM per reaction and mix with P2
	RdRp_SARSr-R	CARATGTTAAASACACTATTAGCATA	Use 800 nM per reaction
E gene	E_Sarbeco_F	ACAGGTACGTTAATAGTTAATAGCGT	Use 400 nm per reaction
	E_Sarbeco_P1	FAM-ACACTAGCCATCCTTACTGCGCTTCG-BBQ	Use 200 nm per reaction
	E_Sarbeco_R	ATATTGCAGCAGTACGCACACA	Use 400 nm per reaction
N gene	N_Sarbeco_F	CACATTGGCACCCGCAATC	Use 600 nm per reaction
	N_Sarbeco_P	FAM-ACTTCCTCAAGGAACAACATTGCCA-BBQ	Use 200 nm per reaction
	N_Sarbeco_R	GAGGAACGAGAAGAGGCTTG	Use 800 nm per reaction

W is A/T; R is G/A; M is A/C; S is G/C. FAM: 6-carboxyfluorescein; BBQ: blackberry quencher.
^b Optimised concentrations are given in nanomol per litre (nM) based on the final reaction mix, e.g. 1.5 µL of a 10 µM primer stock solution per 25 µL total reaction volume yields a final concentration of 600 nM as indicated in the table.

The WHO-protocol (Figure 1), which directly derives from the Corman-Drosten paper, concludes that in order to confirm the presence of SARS-CoV-2, two control genes (the E- and the RdRp-genes) must be identified in the assay. It should be noted, that the RdPd-gene has one uncertain position (“wobbly”) in the forward-primer (R=G/A), two uncertain positions in the reverse-primer (R=G/A; S=G/C) and it has three uncertain positions in the RdRp-probe (W=A/T; R=G/A; M=A/C). So, two different forward primers, four different reverse primers, and eight distinct probes can be synthesized for the RdPd-gene. Together, there are 64 possible combinations of primers and probes!

The Corman-Drosten paper further identifies a third gene which, according to the WHO protocol, was not further validated and deemed unnecessary: *“Of note, the N gene assay also performed well but was not subjected to intensive further validation because it was slightly less sensitive.”*

This was an unfortunate omission as it would be best to use all three gene PCRs as

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

confirmatory assays, and this would have resulted in an almost sufficient virus RNA detection diagnostic tool protocol. Three confirmatory assay-steps would at least minimize-out errors & uncertainties at every fold-step in regards to “Wobbly”-spots. (Nonetheless, the protocol would still fall short of any “good laboratory practice”, when factoring in all the other design-errors).

As it stands, the N gene assay is regrettably neither proposed in the WHO-recommendation (Figure 1) as a mandatory and crucial third confirmatory step, nor is it emphasized in the Corman-Drosten paper as important optional reassurance “for a routine workflow” (Table 2).

Consequently, in nearly all test procedures worldwide, merely 2 primer-matches were used instead of all three. This oversight renders the entire test-protocol useless with regards to delivering accurate test-results of real significance in an ongoing pandemic.

Background

We used known SARS- and SARS-related coronaviruses (bat viruses from our own studies as well as literature sources) to generate a non-redundant alignment (excerpts shown in Annex). We designed candidate diagnostic RT-PCR assays before release of the first sequence of 2019-nCoV. Upon sequence release, the following assays were selected based on their matching to 2019-nCoV as per inspection of the sequence alignment and initial evaluation (Figures 1 and 2).

All assays can use SARS-CoV genomic RNA as positive control. Synthetic control RNA for 2019-nCoV E gene assay is available via EVAg. Synthetic control for 2019-nCoV RdRp is expected to be available via EVAg from Jan 21st onward.

First line screening assay: E gene assay

Confirmatory assay: RdRp gene assay

Figure 1: The N-Gene confirmatory-assay is neither emphasized as necessary third step in the official WHO Drosten-Corman protocol-recommendation [8] nor is it required as a crucial step for higher test-accuracy in the Eurosurveillance publication.

1c) Erroneous GC-content (discussed in 2c, together with annealing temperature (T_m))

1d) Detection of viral genes

RT-PCR is not recommended for primary diagnostics of infection. This is why the RT-PCR Test

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

used in clinical routine for detection of COVID-19 is not indicated for COVID-19 diagnosis on a regulatory basis.

“Clinicians need to recognize the enhanced accuracy and speed of the molecular diagnostic techniques for the diagnosis of infections, but also to understand their limitations. Laboratory results should always be interpreted in the context of the clinical presentation of the patient, and appropriate site, quality, and timing of specimen collection are required for reliable test results”. [9]

However, it may be used to help the physician’s differential diagnosis when he or she has to discriminate between different infections of the lung (Flu, Covid-19 and SARS have very similar symptoms). For a confirmative diagnosis of a specific virus, at least 3 specific primer pairs must be applied to detect 3 virus-specific genes. Preferably, these target genes should be located with the greatest distance possible in the viral genome (opposite ends included). Although the Corman-Drosten paper describes 3 primers, these primers only cover roughly half of the virus’ genome. This is another factor that decreases specificity for detection of intact COVID-19 virus RNA and increases the quote of false positive test results.

Therefore, even if we obtain three positive signals (i.e. the three primer pairs give 3 different amplification products) in a sample, this does not prove the presence of a virus. A better primer design would have terminal primers on both ends of the viral genome. This is because the whole viral genome would be covered and three positive signals can better discriminate between a complete (and thus potentially infectious) virus and fragmented viral genomes (without infectious potency). In order to infer anything of significance about the infectivity of the virus, the Orf1 gene, which encodes the essential replicase enzyme of SARS-CoV-1 and SARS-CoV-2 viruses, should have been included as a target (Figure 2). The positioning of the targets in the region of the viral genome that is most heavily and variably transcribed is another weakness of the protocol.

Kim *et al.* demonstrate a highly variable 3’ expression of subgenomic RNA in Sars-CoV-2 [23]. These RNAs are actively monitored as signatures for asymptomatic and non-infectious patients [10]. It is highly questionable to screen a population of asymptomatic people with qPCR primers that have 6 base pairs primer-dimer on the 3 prime end of a primer (Figure 3). Apparently the WHO recommends these primers. We tested all the wobble derivatives from

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

the Corman-Drosten paper with ThermoFisher's primer dimer web tool [11]. The RdRp forward primer has 6bp 3prime homology with Sarbeco E Reverse. At high primer concentrations this is enough to create inaccuracies.

Of note: There is a perfect match of one of the N primers to a clinical pathogen (*Pantoea*), found in immuno-compromised patients. The reverse primer hits *Pantoea* as well but not in the same region (Figure 3).

These are severe design errors, since the test cannot discriminate between the whole virus and viral fragments. The test cannot be used as a diagnostic for SARS-CoV-2 viruses.

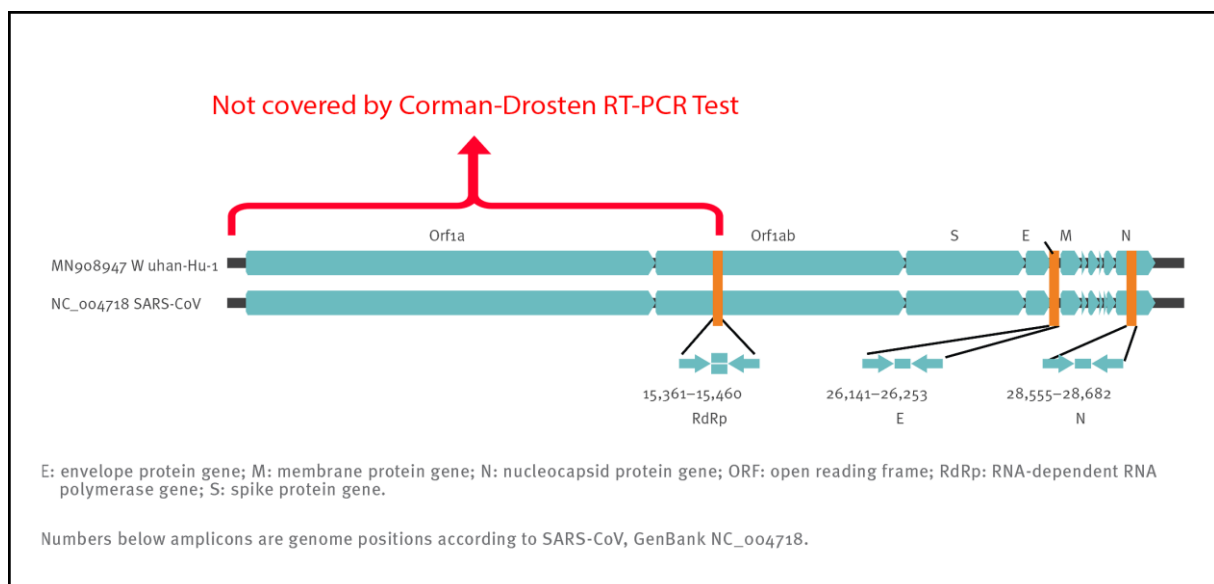


Figure 2: Relative positions of amplicon targets on the SARS-CoV-1 coronavirus and the 2019 novel coronavirus genome. ORF: open reading frame; RdRp: RNA-dependent RNA polymerase. Numbers below amplicon are genome positions according to SARS-CoV-1, NC_004718 [1];

Cross Primer Dimers:

Corman_RdRp_SARs_F1 with Corman_E_Sarbeco_R
Corman_RdRp_SARs_F1
5-gtgaatggtcatgtgtggcgg->
|||||
<-acacacgcatgacgacgttata-5

Corman_RdRp_SARs_F2 with Corman_E_Sarbeco_R
Corman_RdRp_SARs_F2
5-gtgagatggtcatgtgtggcgg->
|||||
<-acacacgcatgacgacgttata-5

> **Corman_N_Sarbeco_F**
CACATTGGCACCCGCAATC

Pantoea agglomerans strain ASB05 chromosome, complete genome
Sequence ID: [CP046722.1](#) Length: 4022781 Number of Matches: 2

Range 1: 2326019 to 2326037 [GenBank](#) [Graphics](#) ▼ Next Match

Score	Expect	Identities	Gaps	Strand
38.2 bits(19)	2.2	19/19(100%)	0/19(0%)	Plus/Plus
Query 1		CACATTGGCACCCGCAATC 19		
Sbjct 2326019		CACATTGGCACCCGCAATC 2326037		

Figure 3: A test with Thermofischer’s primer dimer web tool reveals that the RdRp forward primer has a 6bp 3’ prime homology with Sarbeco E Reverse (left box). Another test reveals that there is a perfect match for one of the N-primers to a clinical pathogen (*Pantoea*) found in immuno-compromised patients (right box).

2. Reaction temperatures

2a) DNA melting temperature (>92°).

Adequately addressed in the Corman-Drosten paper.

2b) DNA amplification temperature.

Adequately addressed in the Corman-Drosten paper.

2c) Erroneous GC-contents and Tm

The annealing-temperature determines at which temperature the primer attaches/detaches from the target sequence. For an efficient and specific amplification, GC content of primers should meet a minimum of 40% and a maximum of 60% amplification. As indicated in table 3, three of the primers described in the Corman-Drosten paper are not within the normal range for GC-content. Two primers (RdRp_SARSr_F and RdRp_SARSr_R) have unusual and very low GC-values of 28%-31% for all possible variants of wobble bases, whereas primer E_Sarbeco_F has a GC-value of 34.6% (Table 3 and second panel of Table 3).

It should be noted that the GC-content largely determines the binding to its specific target due to its three hydrogen bonds in base pairing. Thus, the lower the GC-content of the primer, the lower its binding-capability to its specific target gene sequence (i.e. the gene to

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

be detected). This means for a target-sequence to be recognized we have to choose a temperature which is as close as possible to the actual annealing-temperature (best practise-value) for the primer not to detach again, while at the same time specifically selecting the target sequence.

If the T_m -value is very low, as observed for all wobbly-variants of the RdRp reverse primers, the primers can bind non-specifically to several targets, decreasing specificity and increasing potential false positive results.

The annealing temperature (T_m) is a crucial factor for the determination of the specificity /accuracy of the qPCR procedure and essential for evaluating the accuracy of qPCR-protocols. Best-practice recommendation: Both primers (forward and reverse) should have an almost similar value, preferably the identical value.

We used the freely available primer design software Primer-BLAST [12, 25] to evaluate the best-practise values for all primers used in the Corman-Drosten paper (Table 3). We attempted to find a T_m -value of 60° C, while similarly seeking the highest possible GC%-value for all primers. A maximal T_m difference of 2° C within primer pairs was considered acceptable. Testing the primer pairs specified in the Corman-Drosten paper, we observed a difference of 10° C with respect to the annealing temperature T_m for primer pair1 (RdRp_SARSr_F and RdRp_SARSr_R). *This is a very serious error and makes the protocol useless as a specific diagnostic tool.*

Additional testing demonstrated that only the primer pair designed to amplify the N-gene (N_Sarbeco_F and N_Sarbeco_R) reached the adequate standard to operate in a diagnostic test, since it has a sufficient GC-content and the T_m difference between the primers (N_Sarbeco_F and N_Sarbeco_R) is 1.85° C (below the crucial maximum of 2° C difference). Importantly, this is the gene which was neither tested in the virus samples (Table 2) nor emphasized as a confirmatory test. In addition to highly variable melting temperatures and degenerate sequences in these primers, there is another factor impacting specificity of the procedure: the dNTPs (0.4uM) are 2x higher than recommended for a highly specific amplification. There is additional magnesium sulphate added to the reaction as well. This procedure combined with a low annealing temperature can create non-specific amplifications. When additional magnesium is required for qPCR, specificity of the assay should be further scrutinized.

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

The design errors described here are so severe that it is highly unlikely that specific amplification of SARS-CoV-2 genetic material will occur using the protocol of the Corman-Drosten paper.

Table 3: GC-content of the primers and probes (adapted from Corman-Drosten paper; aberrations from optimized GC-contents are highlighted. Second Panel shows a table-listing of all Primer-BLAST best practices values for all primers and probes used in the Corman-Drosten paper by Prof. Dr. Ulrike Kämmerer & her team

Normal ranges for GC%: 40 - 60%; normal ranges for TM: 55-65°; Best-practise for qPCR in our case: 60° for both primers (reverse & forward)

Assay/use	Oligonucleotide	Sequence*	Concentration*
RdRP gene	RdRp_SARSr-F	GTGARATGGTCATGTGGCGG	Use 600 nM per reaction
	RdRp_SARSr-P2	FAM-CAGGTGGAACCTCATCAGGAGATGC-BBQ	Specific for 2019-nCoV, will not detect SARS-CoV.
	RdRp_SARSr-P1	FAM-CCAGGTGGWACRTCATCMGGTGATGC-BBQ	Use 100 nM per reaction and mix with P1 Pan Sarbeco-Probe will detect 2019-nCoV, SARS-CoV and bat-SARS-related CoVs.
E gene	RdRp_SARSr-R	CARATGTTAAASACACTATTAGCATA	Use 100 nM per reaction and mix with P2
	E_Sarbeco_F	ACAGGTACGTTAATAGTTAATAGCGT	Use 800 nM per reaction
	E_Sarbeco_P1	FAM-ACACTAGCCATCCTTACTGCGCTTCG-BBQ	Use 400 nm per reaction
N gene	E_Sarbeco_R	ATATTGCAGCAGTACGCACACA	Use 200 nm per reaction
	N_Sarbeco_F	CACATTGGCACCCGCAATC	Use 400 nm per reaction
	N_Sarbeco_P	FAM-ACTTCTCTCAAGGAACAACATTGCCA-BBQ	Use 600 nm per reaction
	N_Sarbeco_R	GAGGAACGAGAAGAGGCTTG	Use 200 nm per reaction
			Use 800 nm per reaction

* W is A/T; R is G/A; M is A/C; S is G/C. FAM: 6-carboxyfluorescein; BBQ: blackberry quencher.
* Optimised concentrations are given in nanomol per litre (nM) based on the final reaction mix, e.g. 3.5 µL of a 10 µM primer stock solution per 25 µL total reaction volume yields a final concentration of 600 nM as indicated in the table.

Primer pairs	Sequence (5'-3')	GC Template strand	TM Length	Start	Stop	Tm	GC%	Self 5' complementarity	Self 3' complementarity	Product length (bp)
E_Sarbeco_F	ACAGGTACGTTAATAGTTAATAGCGT	Plus	26	26269	26294	58.29	34.62	8.00	8.00	113
E_Sarbeco_R	ATATTGCAGCAGTACGCACACA	Minus	22	26381	26360	60.93	45.45	7.00	1.00	
N-Sarbeco_F	CACATTGGCACCCGCAATC	Plus	19	28706	28724	60.15	57.89	4.00	0.00	128
N-Sarbeco_R	GAGGAACGAGAAGAGGCTTG	Minus	20	28833	28814	58.00	55.00	3.00	1.00	
RdRp_SARSr-F	GTGARATGGTCATGTGGCGG		22			63.74	59.09	4.00		
RdRp_SARSr-R	CARATGTTAAASACACTATTAGCATA		25			53.56	28.00	7.00		to be added in next version
If R= G and S= G	GTGAGATGGTCATGTGGCGG CAGATGTTAAAGACACTATTAGCATA		22 26			63.74 55.22	59.09 30.77	4.00 7.00	1.00 5.00	
If R= G and S= C	GTGAGATGGTCATGTGGCGG CAGATGTTAAACACACTATTAGCATA		22 26			63.74 55.68	59.09 30.77	4.00 7.00	1.00 2.00	not found in the Sequence
If R= A and S= G	GTGAAATGGTCATGTGGCGG CAAATGTTAAAGACACTATTAGCATA		22 26			62.58 54.23	54.55 26.92	4.00 7.00	1.00 5.00	
If R= A and S= C	GTGAAATGGTCATGTGGCGG CAAATGTTAAACACACTATTAGCATA		22 26			62.58 54.69	54.55 26.92	4.00 7.00	1.00 2.00	
Probes:										
RdRp-SARSr-P2	CAGGTGGAACCTCATCAGGAGATGC		25			64.89	56.00	6.00	5.00	
RdRp-SARSr-P1	CCAGGTGGWACRTCATCMGGTGATGC									
E-Sarbeco-P1	ACACTAGCCATCCTTACTGCGCTTCG		26			66.78	53.85	4.00	2.00	
N-Sarbeco-P	ACTTCTCAAGGAACAACATTGCCA		25			63.15	44.00	8.00	3.00	

3. The number of amplification cycles

It should be noted that there is no mention anywhere in the Corman-Drosten paper of a test being positive or negative, or indeed what defines a positive or negative result. These types of virological diagnostic tests must be based on a SOP, including a validated and fixed number of PCR cycles (Ct value) after which a sample is deemed positive or negative. The maximum reasonably reliable Ct value is 30 cycles. Above a Ct of 35 cycles, rapidly increasing numbers of false positives must be expected .

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

PCR data evaluated as positive after a Ct value of 35 cycles are completely unreliable.

Citing Jaafar *et al.* 2020 [3]: “At Ct = 35, the value we used to report a positive result for PCR, <3% of cultures are positive.” In other words, there was no successful virus isolation of SARS-CoV-2 at those high Ct values.

Further, scientific studies show that only non-infectious (dead) viruses are detected with Ct values of 35 [22].

Between 30 and 35 there is a grey area, where a positive test cannot be established with certainty. This area should be excluded. Of course, one could perform 45 PCR cycles, as recommended in the Corman-Drosten WHO-protocol (Figure 4), but then you also have to define a reasonable Ct-value (which should not exceed 30). But an analytical result with a Ct value of 45 is scientifically and diagnostically absolutely meaningless (a reasonable Ct-value should not exceed 30). All this should be communicated very clearly. It is a significant mistake that the Corman-Drosten paper does not mention the maximum Ct value at which a sample can be unambiguously considered as a positive or a negative test-result. This important cycle threshold limit is also not specified in any follow-up submissions to date.

3. Discriminatory assay		
RdRp assay:		
MasterMix:	Per reaction	
H ₂ O (RNase free)	1.1 µl	
2x Reaction mix*	12.5 µl	
MgSO ₄ (50mM)	0.4 µl	
BSA (1 mg/ml)**	1 µl	
Primer RdRP_SARSr-F2 (10 µM stock solution)	1.5 µl	GTGARATGGTCATGTGTGGCCG
Primer RdRP_SARSr-R1 (10 µM stock solution)	2 µl	CARATGTTAAASACACTATTAGCATA
Probe RdRP_SARSr-P2 (10 µM stock solution)	0.5 µl	FAM-CAGGTGGAACCTCATCAGGAGATGC-BBQ
SSIII/Taq EnzymeMix*	1 µl	
Total reaction mix	20 µl	
Template RNA, add	5 µl	
Total volume	25 µl	

* Thermo Fischer/Invitrogen: SuperScriptIII OneStep RT-PCR System with Platinum® Taq DNA Polymerase
** MgSO₄ (50 mM) [Sigma], This component is not provided with the OneStep RT-PCR kit
*** non-acetylated [Roche].

Cycler:
55°C 10'
94°C 3'
94°C 15"
58°C 30" 45x

Figure 4: RT-PCR Kit recommendation in the official Corman-Drosten WHO-protocol [8]. Only a “Cycler”-value (cycles) is to be found without corresponding and scientifically reasonable Ct (Cutoff-value). This or any other cycles-value is nowhere to be found in the actual Corman-Drosten paper.

4. Biomolecular validations

To determine whether the amplified products are indeed SARS-CoV-2 genes, biomolecular validation of amplified PCR products is essential. For a diagnostic test, this validation is an absolute must.

Validation of PCR products should be performed by either running the PCR product in a 1% agarose-EtBr gel together with a size indicator (DNA ruler or DNA ladder) so that the size of the product can be estimated. The size must correspond to the calculated size of the amplification product. But it is even better to sequence the amplification product. The latter will give 100% certainty about the identity of the amplification product. Without molecular validation one can not be sure about the identity of the amplified PCR products. Considering the severe design errors described earlier, the amplified PCR products can be anything.

Also not mentioned in the Corman-Drosten paper is the case of small fragments of qPCR (around 100bp): It could be either 1,5% agarose gel or even an acrylamide gel.

The fact that these PCR products have not been validated at molecular level is another striking error of the protocol, making any test based upon it useless as a specific diagnostic tool to identify the SARS-CoV-2 virus.

5. Positive and negative controls to confirm/refute specific virus detection.

The unconfirmed assumption described in the Corman-Drosten paper is that SARS-CoV-2 is the only virus from the SARS-like beta-coronavirus group that currently causes infections in humans. The sequences on which their PCR method is based are *in silico* sequences, supplied by a laboratory in China [23], because at the time of development of the PCR test no control material of infectious (“live”) or inactivated SARS-CoV-2 was available to the authors. The PCR test was therefore designed using the sequence of the known SARS-CoV-1 as a control material for the Sarbeco component (Dr. Meijer, co-author Corman-Drosten paper in an email exchange with Dr. Peter Borger) [2].

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

All individuals testing positive with the RT-PCR test, as described in the Corman-Drosten paper, are assumed to be positive for SARS-CoV-2 infections. There are three severe flaws in their assumption. First, a positive test for the RNA molecules described in the Corman-Drosten paper cannot be equated to “infection with a virus”. A positive RT-PCR test merely indicates the presence of viral RNA molecules. As demonstrated under point 1d (above), the Corman-Drosten test was not designed to detect the full-length virus, but only a fragment of the virus. We already concluded that this classifies the test as unsuitable as a diagnostic test for SARS-virus infections.

Secondly and of major relevance, the functionality of the published RT-PCR Test was not demonstrated with the use of a positive control (isolated SARS-CoV-2 RNA) which is an essential scientific gold standard.

Third, the Corman-Drosten paper states:

“To show that the assays can detect other bat-associated SARS-related viruses, we used the E gene assay to test six bat-derived faecal samples available from Drexler et al. [...] und Muth et al. [...]. These virus-positive samples stemmed from European rhinolophid bats. Detection of these phylogenetic outliers within the SARS-related CoV clade suggests that all Asian viruses are likely to be detected. This would, theoretically, ensure broad sensitivity even in case of multiple independent acquisitions of variant viruses from an animal reservoir.”

This statement demonstrates that the E gene used in RT-PCR test, as described in the Corman-Drosten paper, is not specific to SARS-CoV-2. The E gene primers also detect a broad spectrum of other SARS viruses.

The genome of the coronavirus is the largest of all RNA viruses that infect humans and they all have a very similar molecular structure. Still, SARS-CoV-1 and SARS-CoV-2 have two highly specific genetic fingerprints, which set them apart from the other coronaviruses. First, a unique fingerprint-sequence (KTFPPTEPKDKKKK) is present in the N-protein of SARS-CoV-1 and SARS-CoV-2 [13,14,15]. Second, both SARS-CoV-1 and SARS-CoV-2 do not contain the HE protein, whereas all other coronaviruses possess this gene [13, 14]. So, in order to specifically detect a SARS-CoV-1 and SARS-CoV-2 PCR product the above region in the N gene should have been chosen as the amplification target. A reliable diagnostic test should focus on this specific region in the N gene as a confirmatory test. The PCR for this N gene was not further validated nor recommended as a test gene by the Drosten-Corman paper, because of

being “not so sensitive” with the SARS-CoV original probe [1].

Furthermore, the absence of the HE gene in both SARS-CoV-1 and SARS-CoV-2 makes this gene the ideal negative control to exclude other coronaviruses. The Corman-Drosten paper does not contain this negative control, nor does it contain any other negative controls. The PCR test in the Corman-Drosten paper therefore contains neither a unique positive control nor a negative control to exclude the presence of other coronaviruses. This is another major design flaw which classifies the test as unsuitable for diagnosis.

6. Standard Operational Procedure (SOP) is not available

There should be a Standard Operational Procedure (SOP) available, which unequivocally specifies the above parameters, so that all laboratories are able to set up the identical same test conditions. To have a validated universal SOP is essential, because it facilitates data comparison within and between countries. It is very important to specify all primer parameters unequivocally. We note that this has not been done. Further, the Ct value to indicate when a sample should be considered positive or negative is not specified. It is also not specified when a sample is considered infected with SARS-CoV viruses. As shown above, the test cannot discern between virus and virus fragments, so the Ct value indicating positivity is crucially important. This Ct value should have been specified in the Standard Operational Procedure (SOP) and put on-line so that all laboratories carrying out this test have exactly the same boundary conditions. It points to flawed science that such an SOP does not exist. The laboratories are thus free to conduct the test as they consider appropriate, resulting in an enormous amount of variation. Laboratories all over Europe are left with a multitude of questions; which primers to order? which nucleotides to fill in the undefined places? which Tm value to choose? How many PCR cycles to run? At what Ct value is the sample positive? And when is it negative? And how many genes to test? Should all genes be tested, or just the E and RpRd gene as shown in Table 2 of the Corman-Drosten paper? Should the N gene be tested as well? And what is their negative control? What is their positive control? The protocol as described is unfortunately very vague and erroneous in its design that one can go in dozens of different directions. There does not appear to be any standardization nor an SOP, so it is not clear how this test can be implemented.

7. Consequences of the errors described under 1-5: false positive results.

The RT-PCR test described in the Corman-Drosten paper contains so many molecular biological design errors (see 1-5) that it is not possible to obtain unambiguous results. It is inevitable that this test will generate a tremendous number of so-called “false positives”. The definition of false positives is a negative sample, which initially scores positive, but which is negative after retesting with the same test. False positives are erroneous positive test-results, i.e. negative samples that test positive. And this is indeed what is found in the Corman-Drosten paper. On page 6 of the manuscript PDF the authors demonstrate, that even under well-controlled laboratory conditions, a considerable percentage of false positives is generated with this test:

“In four individual test reactions, weak initial reactivity was seen however they were negative upon retesting with the same assay. These signals were not associated with any particular virus, and for each virus with which initial positive reactivity occurred, there were other samples that contained the same virus at a higher concentration but did not test positive. Given the results from the extensive technical qualification described above, it was concluded that this initial reactivity was not due to chemical instability of real-time PCR probes and most probably to handling issues caused by the rapid introduction of new diagnostic tests and controls during this evaluation study.” [1]

The first sentence of this excerpt is clear evidence that the PCR test described in the Corman-Drosten paper generates false positives. Even under the well-controlled conditions of the state-of-the-art Charité-laboratory, 4 out of 310 primary-tests are false positives per definition. Four negative samples initially tested positive, then were negative upon retesting. This is the classical example of a false positive. In this case the authors do not identify them as false positives, which is intellectually dishonest.

Another telltale observation in the excerpt above is that the authors explain the false positives away as "handling issues caused by the rapid introduction of new diagnostic tests". Imagine the laboratories that have to introduce the test without all the necessary information normally described in an SOP.

8. The Corman-Drosten paper was not peer-reviewed

Before formal publication in a scholarly journal, scientific and medical articles are traditionally certified by “peer review.” In this process, the journal’s editors take advice from various experts (“referees”) who have assessed the paper and may identify weaknesses in its assumptions, methods, and conclusions. Typically a journal will only publish an article once the editors are satisfied that the authors have addressed referees’ concerns and that the data presented supports the conclusions drawn in the paper.” This process is as well described for Eurosurveillance [16].

The Corman-Drosten paper was submitted to Eurosurveillance on January 21st 2020 and accepted for publication on January 22nd 2020. On January 23rd 2020 the paper was online. On January 13th 2020 version 1-0 of the protocol was published at the official WHO website [17], updated on January 17th 2020 as document version 2-1 [18], even before the Corman-Drosten paper was published on January 23rd at Eurosurveillance.

Normally, peer review is a time-consuming process since at least two experts from the field have to critically read and comment on the submitted paper. In our opinion, this paper was not peer-reviewed. Twenty-four hours are simply not enough to carry out a thorough peer review. Our conclusion is supported by the fact that a tremendous number of very serious design flaws were found by us, which make the PCR test completely unsuitable as a diagnostic tool to identify the SARS-CoV-2 virus. Any molecular biologist familiar with RT-PCR design would have easily observed the grave errors present in the Corman-Drosten paper before the actual review process. We asked Eurosurveillance on October 26th 2020 to send us a copy of the peer review report. To date, we have not received this report and in a letter dated November 18th 2020, the ECDC as host for Eurosurveillance declined to provide access without providing substantial scientific reasons for their decision. On the contrary, they write that “disclosure would undermine the purpose of scientific investigations.” [24].

9. Authors as the editors

A final point is one of major concern. It turns out that two authors of the Corman-Drosten paper, Christian Drosten and Chantal Reusken, are also members of the editorial board of this journal [19]. Hence there is a severe conflict of interest which strengthens suspicions

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

that the paper was not peer-reviewed. It has the appearance that the rapid publication was possible simply because the authors were also part of the editorial board at Eurosurveillance. This practice is categorized as compromising scientific integrity .

SUMMARY CATALOGUE OF ERRORS FOUND IN THE PAPER

The Corman-Drosten paper contains the following specific errors:

1. There exists no specified reason to use these extremely high concentrations of primers in this protocol. The described concentrations lead to increased nonspecific bindings and PCR product amplifications, making the test unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus.
2. Six unspecified wobbly positions will introduce an enormous variability in the real world laboratory implementations of this test; the confusing nonspecific description in the Corman-Drosten paper is not suitable as a Standard Operational Protocol making the test unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus.
3. The test cannot discriminate between the whole virus and viral fragments. Therefore, the test cannot be used as a diagnostic for intact (infectious) viruses, making the test unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus and make inferences about the presence of an infection.
4. A difference of 10° C with respect to the annealing temperature T_m for primer pair1 (RdRp_SARSr_F and RdRp_SARSr_R) also makes the test unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus.
5. A severe error is the omission of a Ct value at which a sample is considered positive and negative. This Ct value is also not found in follow-up submissions making the test unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus.

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

6. The PCR products have not been validated at the molecular level. This fact makes the protocol useless as a specific diagnostic tool to identify the SARS-CoV-2 virus.
7. The PCR test contains neither a unique positive control to evaluate its specificity for SARS-CoV-2 nor a negative control to exclude the presence of other coronaviruses, making the test unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus.
8. The test design in the Corman-Drosten paper is so vague and flawed that one can go in dozens of different directions; nothing is standardized and there is no SOP. This highly questions the scientific validity of the test and makes it unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus.
9. Most likely, the Corman-Drosten paper was not peer-reviewed making the test unsuitable as a specific diagnostic tool to identify the SARS-CoV-2 virus.
10. We find severe conflicts of interest for at least four authors, in addition to the fact that two of the authors of the Corman-Drosten paper (Christian Drosten and Chantal Reusken) are members of the editorial board of Eurosurveillance. A conflict of interest was added on July 29 2020 (Olfert Landt is CEO of TIB-Molbiol; Marco Kaiser is senior researcher at GenExpress and serves as scientific advisor for TIB-Molbiol), that was not declared in the original version (and still is missing in the PubMed version); TIB-Molbiol is the company which was “the first” to produce PCR kits (Light Mix) based on the protocol published in the Corman-Drosten manuscript, and according to their own words, they distributed these PCR-test kits before the publication was even submitted [20]; further, Victor Corman & Christian Drosten failed to mention their second affiliation: the commercial test laboratory “Labor Berlin”. Both are responsible for the virus diagnostics there [21] and the company operates in the realm of real time PCR-testing.

CONCLUSION

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

In light of our re-examination of the test protocol to identify SARS-CoV-2 described in the Corman-Drosten paper we have identified concerning errors and inherent fallacies which render the SARS-CoV-2 PCR test useless.

The decision as to which test protocols are published and made widely available lies squarely in the hands of Eurosurveillance. A decision to recognise the errors apparent in the Corman-Drosten paper has the benefit to greatly minimise human cost and suffering going forward. Is it not in the best interest of Eurosurveillance to retract this paper? Our conclusion is clear. In the face of all the tremendous PCR-protocol design flaws and errors described here, we conclude: There is not much of a choice left in the framework of scientific integrity and responsibility.

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

References

[1] Corman Victor M, Landt Olfert, Kaiser Marco, Molenkamp Richard, Meijer Adam, Chu Daniel KW, Bleicker Tobias, Brünink Sebastian, Schneider Julia, Schmidt Marie Luisa, Mulders Daphne GJC, Haagmans Bart L, van der Veer Bas, van den Brink Sharon, Wijsman Lisa, Goderski Gabriel, Romette Jean-Louis, Ellis Joanna, Zambon Maria, Peiris Malik, Goossens Herman, Reusken Chantal, Koopmans Marion PG, Drosten Christian. Detection of 2019 novel coronavirus (2019-nCoV) by real-time RT-PCR. Euro Surveill. 2020;25(3):pii=2000045.

<https://doi.org/10.2807/1560-7917.ES.2020.25.3.2000045>

[2] Email communication between Dr. Peter Borger & Dr. Adam Meijer: Supplementary Material

[3] Jafaar *et al.*, Correlation Between 3790 Quantitative Polymerase Chain Reaction–Positives Samples and Positive Cell Cultures, Including 1941 Severe Acute Respiratory Syndrome Coronavirus 2 Isolates

<https://academic.oup.com/cid/advance-article/doi/10.1093/cid/ciaa1491/5912603>

[4] BBC, January 21st 2020: <https://www.bbc.com/news/world-asia-china-51185836>; archive: <https://archive.is/0qRmZ>

[5] Google Analytics - COVID19-deaths worldwide: <https://bit.ly/3fndemJ>; archive: <https://archive.is/PpgEE>

[6] Laboratory testing for COVID-19 Emergency Response Technical Centre, NIVD under China CDC March 15th, 2020:

<http://www.chinacdc.cn/en/COVID19/202003/P020200323390321297894.pdf>

[7] Real-Time PCR Handbook Life Technologies

(<https://www.thermofisher.com/content/dam/LifeTech/global/Forms/PDF/real-time-pcr-handbook.pdf>)

Nolan T, Huggett J, Sanchez E. Good practice guide for the application of quantitative PCR (qPCR) First Edition 2013

[8] Trestan Pillonel *et al.*, Letter to the editor: SARS-CoV-2 detection by real-time RT-PCR:

<https://www.ncbi.nlm.nih.gov/pmc/articles/PMC7268274/>

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

[9] Kurkela, Satu, and David WG Brown. "Molecular-diagnostic techniques." *Medicine* 38.10 (2009): 535-540.

[10] Wolfel *et al.*, Virological assessment of hospitalized patients with COVID-2019

<https://www.nature.com/articles/s41586-020-2196-x>

[11] Thermofischer Primer Dimer Web Tool:

<https://www.thermofisher.com/us/en/home/brands/thermo-scientific/molecular-biology/molecular-biology-learning-center/molecular-biology-resource-library/thermo-scientific-web-tools/multiple-primer-analyzer.html>

[12] Primer-BLAST, NCBI - National Center for Biotechnology Information:

<https://www.ncbi.nlm.nih.gov/tools/primer-blast/>

[13] Marra MA, Steven JM, Caroline RA, Robert AH, Angela BW *et al.* (2003) *Science*. The Genome sequence of the SARS-associated coronavirus. *Science* 300(5624): 1399-1404.

[14] Severe acute respiratory syndrome coronavirus 2 isolate Wuhan-Hu-1, complete genome: <https://www.ncbi.nlm.nih.gov/nuccore/MN908947>

[15] Borger P. A SARS-like Coronavirus was expected but nothing was done to be prepared.

Am J Biomed Sci Res 2020. <https://biomedgrid.com/pdf/AJBSR.MS.ID.001312.pdf>

[https://www.researchgate.net/publication/341120750_A_SARS-](https://www.researchgate.net/publication/341120750_A_SARS-like_Coronavirus_was_Expected_but_nothing_was_done_to_be_Prepared)

[like Coronavirus was Expected but nothing was done to be Prepared](https://www.researchgate.net/publication/341120750_A_SARS-like_Coronavirus_was_Expected_but_nothing_was_done_to_be_Prepared); archive:

<https://archive.is/i76Hu>

[16] Eurosurveillance paper evaluation / review process:

<https://www.eurosurveillance.org/evaluation>

[17] Official recommendation of the Corman-Drosten protocol & manuscript by the WHO, published on January 13th 2020 as version 1.0 of the document:

[https://www.who.int/docs/default-source/coronaviruse/wuhan-virus-assay-](https://www.who.int/docs/default-source/coronaviruse/wuhan-virus-assay-v1991527e5122341d99287a1b17c111902.pdf)

[v1991527e5122341d99287a1b17c111902.pdf](https://www.who.int/docs/default-source/coronaviruse/wuhan-virus-assay-v1991527e5122341d99287a1b17c111902.pdf); archive: <https://bit.ly/3m3jXVH>

[18] Official WHO-recommendation for the Corman / Drosten RT-qPCR-protocol, which

directly derives from the Eurosurveillance-publication, document-version 2-1, published on

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

17th January 2020: https://www.who.int/docs/default-source/coronaviruse/protocol-v2-1.pdf?sfvrsn=a9ef618c_2

[19] Eurosurveillance Editorial Board, 2020: <https://www.eurosurveillance.org/upload/site-assets/imgs/2020-09-Editorial%20Board%20PDF.pdf>; archive: <https://bit.ly/2TqXBjX>

[20] Instructions For Use LightMix SarbecoV E-gene plus EAV Control, TIB-Molbiol & Roche Molecular Solutions, January 11th 2020:

[https://www.roche-as.es/lm_pdf/MDx_40-0776_96_Sarbeco-E-gene_V200204_09164154001\(1\).pdf](https://www.roche-as.es/lm_pdf/MDx_40-0776_96_Sarbeco-E-gene_V200204_09164154001(1).pdf)

Archive, timestamp - January 11th 2020: <https://archive.is/Vulo5>; archive: <https://bit.ly/3fm9bXH>

[21] Christian Drosten & Victor Corman, responsible for viral diagnostics at Labor Berlin: <https://www.laborberlin.com/fachbereiche/virologie/>; archive: archive.is/CDEUG

[22] Tom Jefferson, Elizabeth Spencer, Jon Brassey, Carl Heneghan Viral cultures for COVID-19 infectivity assessment. Systematic review. Systematic review doi:

<https://doi.org/10.1101/2020.08.04.20167932>;

<https://www.medrxiv.org/content/10.1101/2020.08.04.20167932v4>

[23] Kim *et al.*, The Architecture of SARS-CoV-2 Transcriptome:

<https://www.sciencedirect.com/science/article/pii/S0092867420304062>

[24] ECDC reply to Dr. Peter Borger, 18th November 2020: Supplementary Material

[25] Prof. Dr. Ulrike Kämmerer & team, survey & Primer-BLAST table: Supplementary Material

Additional literature:

Description RT-PCR RKI Germany, on page 10 of this link:

https://www.rki.de/DE/Content/Gesundheitsmonitoring/Gesundheitsberichterstattung/GBEDownloads/JoHM_S5_2020_Studienprotokoll_CORONA_MONITORING_lokal.pdf?blob=publicationFile

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

Author's Affiliations:

- 1) **Dr. Pieter Borger** (MSc, PhD), Molecular Genetics, W+W Research Associate, Lörrach, Germany,
- 2) Rajesh Kumar Malhotra (Artist Alias: **Bobby Rajesh Malhotra**), Former 3D Artist / Scientific Visualizations at CeMM - Center for Molecular Medicine of the Austrian Academy of Sciences (2019-2020), University for Applied Arts - Department for Digital Arts Vienna, Austria
- 3) **Dr. Michael Yeadon** BSc(Hons) Biochem Tox U Surrey, PhD Pharmacology U Surrey. Managing Director, Yeadon Consulting Ltd, former Pfizer Chief Scientist, United Kingdom,
- 4) **Dr. Clare Craig** MA, (Cantab) BM, BCh (Oxon), FRCPath, United Kingdom
- 5) **Kevin McKernan**, BS Emory University, Chief Scientific Officer, founder Medical Genomics, engineered the sequencing pipeline at WIBR/MIT for the Human Genome Project, Invented and developed the SOLiD sequencer, awarded patents related to PCR, DNA Isolation and Sequencing, USA
- 6) **Prof. Dr. Klaus Steger**, Department of Urology, Pediatric Urology and Andrology, Molecular Andrology, Biomedical Research Center of the Justus Liebig University, Giessen, Germany
- 7) **Dr. Paul McSheehy** (BSc, PhD), Biochemist & Industry Pharmacologist, Loerrach, Germany
- 8) **Dr. Lidiya Angelova**, MSc in Biology, PhD in Microbiology, Former researcher at the National Institute of Allergy and Infectious Diseases (NIAID), Maryland, USA
- 9) **Dr. Fabio Franchi**, Former Dirigente Medico (M.D) in an Infectious Disease Ward, specialized in "Infectious Diseases" and "Hygiene and Preventive Medicine", Società Scientifica per il Principio di Precauzione (SSPP), Italy
- 10) **Dr. med. Thomas Binder**, Internist and Cardiologist (FMH), Switzerland
- 11) **Prof. Dr. med. Henrik Ullrich**, specialist Diagnostic Radiology, Chief Medical Doctor at the Center for Radiology of Collm Oschatz-Hospital, Germany
- 12) **Prof. Dr. Makoto Ohashi**, Professor emeritus, PhD in Microbiology and Immunology, Tokushima University, Japan
- 13) **Dr. Stefano Scoglio**, B.Sc. Ph.D., Microbiologist, Nutritionist, Italy
- 14) **Dr. Marjolein Doesburg-van Kleffens** (MSc, PhD), specialist in Laboratory Medicine (clinical chemistry), Maasziekenhuis Pantein, Beugen, the Netherlands
- 15) **Dr. Dorothea Gilbert** (MSc, PhD), PhD Environmental Chemistry and Toxicology. DGI Consulting Services, Oslo, Norway

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

- 16) **Dr. Rainer J. Klement**, PhD. Department of Radiation Oncology, Leopoldina Hospital Schweinfurt, Germany
- 17) **Dr. Ruth Schruefer**, PhD, human genetics/ immunology, Munich, Germany,
- 18) **Dra. Berber W. Pieksma**, General Practitioner, The Netherlands
- 19) **Dr. med. Jan Bonte** (GJ), Consultant Neurologist, the Netherlands
- 20) **Dr. Bruno H. Dalle Carbonare** (Molecular biologist), IP specialist, BDC Basel, Switzerland
- 21) **Dr. Kevin P. Corbett**, MSc Nursing (Kings College London) PhD (London South Bank) Social Sciences (Science & Technology Studies) London, England, UK
- 22) **Prof. Dr. Ulrike Kämmerer**, specialist in Virology / Immunology / Human Biology / Cell Biology, University Hospital Würzburg, Germany

Review Report - Corman-Drosten *et al.*, Eurosurveillance 2020

Author's Contributions:

PB: Planned and conducted the analyses and research, conceptualising the manuscript.

RKM: Planned and conducted the research, conceptualising the figures and manuscript.

MY: Conducted the analyses and research.

KMcK: Conducted the analyses and research, conceptualized the manuscript.

KS: Conducted the analyses and research.

PMcS: Proofreading the analyses and research.

LA: Proofreading the analyses and research.

FF: Proofreading the analyses and research.

TB: Proofreading the analyses and research.

HU: Proofreading the analyses and research.

MO: Proofreading the analyses and research.

SS: Proofreading the analyses and research.

MDvK: Proofreading the analyses and research.

DG: Proofreading the analyses and research.

RJK: Proofreading the analyses and research.

RS: Proofreading the analyses and research, and the manuscript.

BWK: Proofreading the analyses and research.

RvV: Proofreading the analyses and research.

JB: Proofreading the analyses and research.

KC: Proofreading the analyses and research.

UK: Planned and conducted the analyses and research, conceptualising the manuscript.

Acknowledgement:

We are grateful to Saji N Hameed (Environmental Informatics, University of Aizu, Tsuruga, Ikki-machi, Aizuwakamatsu-shi, Fukushima, Japan) and Howard R. Steen (MA Chem. Eng. Cantab (1969-'73), Former Research Manager, UK) for proofreading our manuscript.